



ESTD 1979

JOTHI X-RAY, E.C.G. & LABORATORY A/c

DIGITAL X-RAY, COMPUTERISED LAB

ISO 9001 CERTIFIED LAB

Ln. P. CHEMBULINGAM, DMLT.,

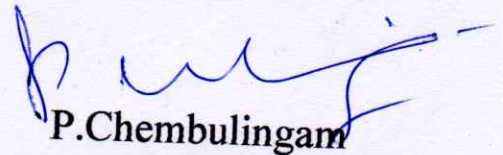
Respected Sir,

You have asked for PT report in NABL portal for our lab. PT report is in process from CMC lab. PT report will get ready only on january month.

Thanking You,

Date : 29-12-2020.

Yours Truly,


P.Chembulingam

LAB MONTHLY SUMMARY

Lab Name **JOTHI X-RAY AND LABORATORY** Lab No **9661**
 Month **October** Year **2020**
 Constituent **Chemistry II**
 Group

Date of Result Entered : 12/10/2020

Date of Report Published : 09/11/2020

Sl.No	Analyte	Method / Principle Name	Analyzer Name	No of Participants	DV	Participants		Your Value	SDI	U
						CV	SD			
1	GLUCOSE I	GOD-POD	Biosystems	232	170.18	9.30	15.83	181 mg/dl	0.68	2.08
2	UREA I	UREASE UV / GLDH	Biosystems	147	108.58	12.42	13.48	99 mg/dl	-0.71	2.22
3	CREATININE I	JAFFE RATE-BLANKED & COMPENSATED/Alkaline Picrate	Biosystems	209	4.71	11.56	0.55	5.5 mg/dl	1.45	0.08
4	T.BILIRUBIN I	EVELYN-MALLOY	Any Analyser	74	3.54	15.84	0.56	6.2 mg/dl	4.75	0.13
5	T-PROTEIN I	BIURET - colorimetric	Biosystems	181	5.02	9.69	0.49	5.4 g/dl	0.78	0.07
6	ALBUMIN I	BCG - colorimetric	Biosystems	188	3.00	9.17	0.28	2.6 g/dl	-1.45	0.04
7	URIC ACID I	ENZYMATIC / URICASE Colorimetric	Biosystems	153	4.36	13.90	0.61	6.2 mg/dl	3.02	0.10
8	CHOLESTEROL I	CHOD-PAP	Biosystems	206	98.79	10.78	10.65	135 mg/dl	3.4	1.48
9	TRIGLYCERIDE I	GPO-PAP / Enzymatic Colorimetric / End Point	Biosystems	189	189.02	10.67	20.16	192 mg/dl	0.15	2.93

SDI Range	Interpretation
Within -1.0 to +1.0	Excellent.
Between ± 1.0 to ± 2.0	Good.
Between ± 2.0 to ± 3.0	Accept with caution. Warning Signal.
Beyond ± 3.0	Unacceptable performance. Action Signal.

Page 1 of 1

Homogeneity and Stability of the sample is passed.

Data in CMC EQAS reports is confidential

Contact details:

Email: clinqc@cmcvellore.ac.in

Contact Number: 0416-2283102

Pamela Christudoss

Dr. Pamela Christudoss
 CMC EQAS Co-Ordinator
 Christian Medical College, Vellore

***** End of Report *****

LAB MONTHLY SUMMARY

Lab Name **JOTHI X-RAY AND LABORATORY** Lab No **9661**
 Month **October** Year **2020**
 Constituent Group **HbA1c**

Date of Result Entered : 12/10/2020

Date of Report Published : 09/11/2020

Sl.No	Analyte	Method / Principle Name	Analyzer Name	No of Participants	DV	Participants		Your Value	SDI	U
						CV	SD			
1	HbA1c	HPLC-ION EXCHANGE	Biorad (D10 / D100 / Variant Turbo)	461	6.60	6.37	0.42	7 %	0.95	0.04

SDI Range	Interpretation
Within -1.0 to +1.0	Excellent.
Between ±1.0 to ±2.0	Good.
Between ±2.0 to ±3.0	Accept with caution. Warning Signal.
Beyond ±3.0	Unacceptable performance. Action Signal.

Page 1 of 1

Homogeneity and Stability of the sample is passed.

Data in CMC EQAS reports is confidential

Contact details:

Email: clinqc@cmcvellore.ac.in

Contact Number: 0416-2283102

Pamela Christudoss

Dr. Pamela Christudoss
CMC EQAS Co-Ordinator
Christian Medical College, Vellore

***** End of Report *****



ESTD 1979

JOTHI X-RAY, E.C.G. & LABORATORY A/c

DIGITAL X-RAY, COMPUTERISED LAB

ISO 9001 CERTIFIED LAB

Ln. P. CHEMBULINGAM, DMLT.,

Clinical Biochemistry

SL. NO	Name of the PTP	Type of Sample	Specific Test	Test Method/techniques	Date of Issue of PT Result	Is Result satisfactory
01	CMC EQAS Cl. Biochemistry, Vellore	Serum	Renal function test- Creatinine	JAFFE	09-11-2020	Yes

Signature

JOTHI X-RAY, ECG & LA

#1, Sivan Kovil Street,

TUTICORIN - 628002

Ph: 2327688 Cell: 9842111104



JOTHI DIAGNOSTIC CENTRE

CORRECTIVE ACTION REPORT

Form No.	JDC / QAS-BC/230
Issue No.	01
Rev No.	0
Issue Date	12.1.2020

CAR No.:

Date: 12-12-2020

Sample

Process

System

Cust. Complaint

Detail of Non-Conformance

We used Biosystems analyser for S. Bilirubin, uric acid & S. Cholesterol. So the result value was high.

Signature:

Date: 12-12-2020

Root Cause Analysis

We started using ~~Biochemistry~~ ^{Biochemistry} Internal Quality Control sample test.

Signature:

Date: 12-12-2020

Corrective Action

We have changed the brand of the Reagent Kit to get correct results for S. Bilirubin, Uric acid & S. Cholesterol.

Signature:

Date: 12-12-2020

Verification & Effectiveness of Corrective Action

We got counselling/clarification from ^{the accurate} manager of the company. Also, we got ^{an} Internal Quality control report.

Signature:

Date: 12-12-2020

Signature

JOTHI X-RAY, ECG & LA

#1, Sivan Kovil Street,

TUTICORIN - 628002

Ph: 2327688 Cell: 9842111104