

TRANSASIA®

Date: 19/06/2021
Effective Date: 19/06/2021

Certificate of Calibration

Customer Name: Narayan Pathology

Model : Automated Hematology Analyzer Sysmex XP-100

Serial No. : A5525

Calibration Done Date: 19.6.21


Next Calibration Due Date On or Before: 18/06/2022

Lab In-charge: _____

This is to certify that the above-mentioned product has been verified of calibration for CBC 5 parameters (WBC, RBC, HGB, HCT and PLT) according to the standard procedures provided by Sysmex Corporation, Japan.

The reference instruments used for value-assignment are managed by the traceability system in Sysmex Corporation and these are traceable to the International Standards, such as ICSH.

Calibration at site performed by
Engineer Name
Designation
Transasia Bio-Medicals Ltd
Location
SASIA BIO-MEDICALS LTD.

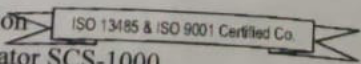incl:
Certificate of Inspection 
Assay Sheet of Calibrator SCS-1000



SASIA BIO-MEDICALS LTD.
HOUSE 8 CHANDIVALLI STUDIO ROAD, MUMBAI-400 072 TEL : 4030 9000 FAX : (022) 4030 9090 / 2857 3030
asia@transasia.co.in
25732223 CHENNAI (044) 28227149 KOLKATA (033) 22157839 BANGALORE (080) 25568044 AHMEDABAD (079) 26407050

Date: 19-06-2021
Effective Date: 19-06-2021


Certificate of Inspection

1. Model: Automated Hematology Analyzer Sysmex XP – 100
2. Serial No.: A5525
3. Calibration Date: 19-06-2021
4. Material used: SCS-1000 (Lot No. 11380525, Expiry date: 20-June-2021)

By comparing your data to the results of the standard counters in Sysmex Corporation, the calibration for CBC 5 parameters using the measurement standard material (SCS-1000) was completed. The calibration result of 5 runs is summarized in the following table. Please refer to the attached sheets for the details.

TRANSASIA BIO-MEDICALS LTD.

transasia.co.in

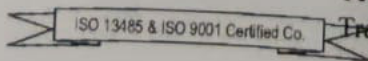

Technical Service Department

Transasia Bio-Medicals Ltd



TRANSASIA BIO-MEDICALS LTD.

TRANSASIA HOUSE, 8 CHANDIVALI STUDIO ROAD, MUMBAI-400 072 TEL : 4030 9000 FAX : (022) 4030 9090 / 2857 3030
transasia@transasia.co.in
(011) 25732223 CHENNAI (044) 28227149 KOLKATA (033) 22157839 BANGALORE (080) 25568044 AHMEDABAD (079) 26407050

TRANSASIA[®]


6. PRECISION STUDY PERFORMED ON THE ANALYZER USING A BLOOD SAMPLE (ORIGINALS ATTACHED)

SMP NO	WBC	RBC	HGB	HCT	PLT
1	8.3	5.10	12.6	40.2	264
2	8.3	5.08	12.7	40.2	265
3	8.3	5.09	12.7	40.2	267
4	8.30	5.14	12.7	40.5	268
5	8.5	5.09	12.7	40.0	260
6	8.3	5.09	12.7	40.0	265
7	8.5	5.15	12.8	40.5	259
8	8.33	5.13	12.8	40.2	239
9	8.30	5.19	12.9	40.6	255
10	8.4	5.17	12.9	40.6	236
Mean	8.35	5.12	12.75	40.30	257.80
SD	0.084	0.039	0.097	0.231	11.419
CV%	1.000	0.753	0.762	0.573	4.430
Acceptable CV%	Within 3.5%	Within 2.0%	Within 1.5%	Within 2.0%	Within 6.0%
Result	PASS	PASS	PASS	PASS	PASS

Technical Service Department
Transasia Bio-Medicals Ltd

TRANSASIA BIO-MEDICALS LTD.

www.transasia.co.in

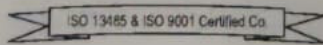

TRANSASIA BIO-MEDICALS LTD.

TRANSASIA HOUSE, 8 CHANDIVALI STUDIO ROAD, MUMBAI-400 072 TEL : 4030 9000 FAX : (022) 4030 9090 / 2857 3030
Email: transasia@transasia.co.in
DELHI (011) 25732223 CHENNAI (044) 28227149 KOLKATA (033) 22157839 BANGALORE (080) 25568044 AHMEDABAD (079) 26407050


TRANSASIA®

5. BACKGROUND CHECK

PARAMETER	RESULT	Range
WBC	0.0	0.3×10^3 /UI or Less
RBC	0.00	0.02×10^6 /uL or Less
HGB	0.0	0.1 g/dL or Less
PLT	0	10×10^3 /uL or Less

IA BIO-MEDICALS LTD.

Technical Service Department
Transasia Bio-Medicals Ltd

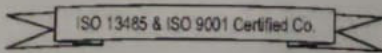

BIO-MEDICALS LTD.

HOUSE, 8 CHANDIVALI STUDIO ROAD, MUMBAI-400 072 TEL : 4030 9000 FAX : (022) 4030 9090 / 2857 3030
@transasia.co.in
32223 CHENNAI (044) 28227149 KOLKATA (033) 22157839 BANGALORE (080) 25568044 AHMEDABAD (079) 26407050

TRANSASIA®

7. CALIBRATION DATA

SMP NO/TIME	WBC	RBC	HGB	HCT	PLT
1	7.10	4.33	11.2	31.00	270
2	7.00	4.32	11.2	31.50	271
3	7.00	4.29	11.2	31.00	268
4	7.20	4.28	11.2	31.06	265
5	7.00	4.34	11.1	31.11	267
MEAN	7.06	4.312	11.18	31.13	268.2
Acceptable Limits	6.69 - 7.29	4.272 - 4.447	11.09 - 11.32	30.67 - 32.05	252.3 - 278.8
Result	PASS	PASS	PASS	PASS	PASS


8. (Traceability System) :

The traceability system of Sysmex Hematology analyzers are shown in attached sheet.

Technical Service Department
Transasia Bio-Medicals Ltd

TRANSASIA BIO-MEDICALS LTD.

www.transasia.co.in

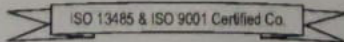

TRANSASIA BIO-MEDICALS LTD.

TRANSASIA HOUSE, 8 CHANDIVALI STUDIO ROAD, MUMBAI-400 072 TEL : 4030 9000 FAX : (022) 4030 9090 / 2857 3030
Email transasia@transasia.co.in
DELHI (011) 25732223 CHENNAI (044) 28227149 KOLKATA (033) 22157839 BANGALORE (080) 25568044 AHMEDABAD (079) 26407050