



POCT SERVICES

Complete Medical Solution

208/207 Transport Nagar, Kanya Board, Madhepura

Transport Nagar, Madhepura, Bihar 800 021, INDIA

Tel: 0947 2412023

Email: poctservices@poctservices.com www.poctservices.com

Date 23 June, 2023

Certificate of Calibration

Name & Address of Customer Pathology Lab Department of JTKMCH Madhepura

City Madhepura State Bihar PIN _____

Phone _____ E-Mail _____

Name of Instrument Selectra PRO M

Type Random Access Fully Automatic Biochemistry Analyser

Serial No: 21-4093

Calibration Date: 21/06/2023

Next Calibration Due 20/6/2024

This is to certify that above said instrument has been validated of hardware calibration for Filters, Aspiration, and Temperature & Lamp according to the procedures provided by Elitech Group Clinical Systems, France


This calibration is carried out by using Standard Operating Procedures (S O P) provided by Elitech Group, shown in the attachment

These instruments conform to CE-IVD & EU directives of use

Calibration carried out on site by - Mr Manish Kumar

Signature & Stamp

Name of Engineer/ Application Specialist - Ajit Mohan Dubey
Asst Sr Manager Application
POCT Services

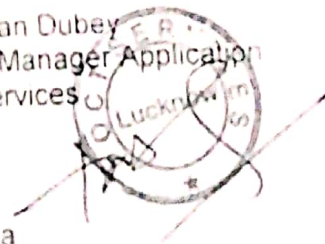

Encls - SOP of Validation/Calibration along with data

[Signature]
7/10/23

प्रभाती चिकित्सा परामर्शकारी
स्वच्छ अस्पताल
एन.ए. रोड, प्रभाती चिकित्सा संस्थान,
मधेपुरा

[Signature]

Nehal



Validation / Calibration - SOP

Selectra ProM

Name of the Customer & Address : Pathology Lab

Address: -----JTKMCH Madhepura

Sr No 21-4093

Status : Under warranty : Validation & Preventive Maintenance

➤ Power Supply

Measure Input power Supply Voltage 229_V (230 V AC \pm 10 V)

Check Earthling: 2.2 V (0 - 5 V)

➤ Ambient temperature: 26 ° C (10 - 35 ° C)

➤ Appearance : Clean (Clean/Dusty)

➤ Bellow Pumps: Open the pump assays and clean it thoroughly.

➤ Analyser Control

Filter: Select the desired position through the Service menu
Filter wheel sets the desired Filter Yes

Filter Status: Needs replacement (Yes/ NO)

ρ 340nm ρ 405nm ρ 505nm ρ 546 nm ρ 578 nm ρ 620 nm ρ 660 nm ρ 700 nm

Note: Filter checked status was ok no need of replacement.

Temperature: Select the desired Options through the Service mode
Temperature OK: Yes

Pump: Select the desired volume through the Service mode
Verify by aspirating the same Quantity : OK

[Signature]
7/10/23

Nikhil



Valve: Select the desired position through the service mode
Valve is energized Yes

Syringes: Check for syringe leakage by physical inspection of syringes
No water leakage Found

Cuvette Drier Block: Check the condition of cuvette drier block by removing the cover of cuvette rotor and lifting the wash arm through service menu. It should be reasonably clean. If dirty please change the drier block

Note:- Condition of cuvette drier block is clean. No need to change

Mixer Belts – Check the elasticity of mixer belts. Should be reasonably good or replace the belts.

Note:- Mixer belts are good no need to change.

Cuvette Rotor Blank : Perform rotor blank and check the OD values of cuvettes. All cuvette blank OD values should be within acceptable range. If required replace the cuvette rotor

Note:- All cuvette blank OD values are in range no need to replace

[Signature]
7/10/23

Nehal



Hardware Calibration of Selectra Pro S/Pro M

➤ Lamp Calibration/Alignment

Lamp Adjustment -

1. Flush the system with distilled water by doing Rotor Blank.
2. Select Adjust Lamp in service menu. Check Value obtained on Display (Adjust the lamp if it is out of 1 800 to 4 000, to as low as possible)

Do not touch lamp !! It may be Hot !!

Lamp alignment Data @ 340 nm wavelength			
Lamp Abs Obtained	Acceptable Range	Alignment	Remarks
3.7229 Abs	1.800 to 4.200abs	Done	Lamp O.D. in acceptable range. No replacement required

➤ Checking the filters

Perform filter check in adjust lamp mode in service menu.

All the arrows must be in Green area. If not, then adjust lamp or replace filter if necessary.

Note :
 When the absorbance value is too low to measure, i.e., the gain is too high in this case, instead of the absorbance value, the value -99999 is shown

Filter (Wavelength)	Gain Range	Gain Achieved	Remarks	Corrective Action
340	0.1 - 3.5	2.7333	OK	Not required
405	0.1 - 2.6	1.5266	OK	Not required
505	0.1 - 2.6	1.0578	OK	Not required
546	0.1 - 2.6	0.5435	OK	Not required
578	0.1 - 2.6	0.7665	OK	Not required
620	0.1 - 1.2	0.6806	OK	Not required
660	0.1 - 0.7	0.6906	OK	Not required
700	0.1 - 0.7	0.6263	OK	Not required
Over all Remarks	Filter gains within acceptable range. No replacement required			

If it is necessary to replace defective filters, please contact service department

[Signature]
 7/10/23

Nehal



➤ **Calibration/Verification of performance of Pipetting system & measuring unit**

- Install dichromate solution on reagent rotor(s) & as sample on sample rotor (Use service disk which has Pre-defined protocol installed for dichromate run)
- Run 10x "Check-S" or 10x "Check-R" as QC samples

Test	Target Value	Target CV [%]	Mean Result	CV [%]
Check-S	0.03(0.060-0.100)	≤2%	0.078	1.378
Check-R	1.75(1.500-2.000)	≤2%	1.573	0.323

Remarks:

Rotor Blank acceptable. Instrument ready for chemical installation & calibration

- Change reagent disk from Service to Standard
- Install the various reagents on reagent rotor(s)
- Install ISE reagents on reagent rotor(s) (If applicable)
- Run Reagent Blanks(s)
- Run Calibrations

Volume calibration of pipettors:-

It is possible to check a predetermined amount of water to check the correct functioning of the pump. Before carrying out this check, the instrument must first carry out a flush routine to ensure that all system tubes are completely filled with water by doing fill system

- 1 Go to Sample syringe full stroke. (For Pro M Model Only)
- 2 Collect the dispensed water. Check the dispensed volume using calibrated pipette (For Pro M only)

Nebel



POCT SERVICES

Complete Medical Solution

298-281, Transport Nagar, Kanpur Road, Adjacent
Transport Nagar Metro Station, Lucknow-226 023 UP, INDIA
Tel: 0522 2433023

E-mail: poc@poctservices.com www.poctservices.com

Pipettor Calibration Data using distilled Water		
Full stroke volume to be dispensed (µl)	Dispensed volume checked and found complying as full stroke volume? (Yes/No)	Remarks
Sample Syringe:		
100	Yes	OK
100	Yes	OK
100	Yes	OK
Reagent syringe:		
1000	Yes	OK
1000	Yes	OK
1000	Yes	OK

Data for volumes other than full stroke:-

This can be verified using pre-determined amount of distilled water in sample/reagent cups and running any dummy program. As soon as the reagent probe/sample probe takes up the sample/reagent, those cups/bottles are taken back and verified for remaining volume using calibrated pipette. Same can be repeated for variable volumes by changing the aspiration volumes in test programmes

Pippettor Calibration Data using distilled Water				
Measured Volume taken in sample cup(µL) (A)	Water to be aspirated by syringe(µL) (B)	Water that should be remaining in cup after aspiration(µL) (C=A-B)	Is the remaining volume inside the cup was found to be the same as in column C? (Yes/No)	Remarks
Reagent syringe:				
5000 µL	300 µL X 3 test = 900 µL	4100 µL	Yes	OK
Sample syringe:				
300 µL	30 µL X 3 test = 90 µL	210 µL	Yes	OK



POCT SERVICES

Complete Medical Solution

298 201, Transport Nagar, Karpur Road, Adjacent
Transport Nagar Metro Station, Lakkur, 726 023, Jh. INDIA

Tel : 0522 2433023

E-mail: poctso@poctservices.com | www.poctservices.com

➤ Temperature Calibration

Select Temperature in Service Menu. It should be $37\text{ }^{\circ}\text{C} \pm 2\text{ }^{\circ}\text{C}$
Verify with temperature Indicator by surface probe in cuvette rotor. If any discrepancy
add the offset of difference in actual & desired temperature.

Temperature Calibration Data				
Displayed Temp	Ref. Range	Temp. Indicator	Temp Offset Required	Temp offset Value
37°C	$37\text{ }^{\circ}\text{C} \pm 2\text{ }^{\circ}\text{C}$	36.7°C	No	0 °C
Remarks	Temp. Calibration OK. No offset required			



➤ Reagent Calibration of the Instrument

Customer is advised to verify the hardware calibration by reagent calibration. Use Elitech Calibrator Elical 2 for the calibration of all parameters.

User can do the same & attach the results in separate sheet with factors after verifying the same with Elitech Elitrol I & Elitrol II controls. All control values should fall within acceptable range. Data sheets of the same should be attached along with this document.

- Switch Off the instrument
- Ensure all the Recommended Spares / Consumables have been replaced. (if not done during PM and required)
- Clean the instrument
- Close the cover

Recommended Spares for replacement : NIL

We hereby certify that Validation have been carried out under the MOU. Hardware Calibration of Lamp, Filters, Temperature & Aspiration (Pump) has been done successfully. Please perform the standardization / Calibration and verify by evaluating controls before processing patient samples.

Next Calibration is due on: 20/06/2024.

Signature of Application Specialist

Place
Date

Nehal