

# Calibration Certificate

**Following voltages have been checked & adjusted on DiaSys SYS 200 Pro (Sl. No. SYS20220401) at Krsnaa Diagnostics – Rajpura, Punjab**

## Calibration Report:

System voltages are in the specified range – Checked by Digital Multimeter are as below:

S.No.	Module	Observed Value	Acceptable Range
1	Power supply board o/p	+24.1 VDC	+24.0 +/- 0.5 VDC
		+05.0 VDC	+05.0 +/- 0.2 VDC
2	Lamp voltage	+12.1 VDC	12.0 +/- 0.5 VDC

System Temperature values are in the specified range –Checked by Digital Thermometer are as below:

S.No.	Module	Observed Value	Acceptable Range
1	Reaction Disk Temp	37.0 C	37.0 C +/- 0.3 C

Photometer Adjustment Values are within range as recommended:

S.No.	Module	Observed Value	Acceptable Range
1	ADC Count for all Wavelengths	10,240 to 10,300	8000 to 18,000

- All other mechanisms have been checked through field service program & Software, found OK.
- Performed calibration for End Point, Kinetic & Fixed Time Kinetic. Run control & sample. Results were satisfactory.
- Next Calibration will be performed on 28/05/2024.

For DiaSys Diagnostics India Pvt. Ltd.



Asmat Sami

Date: 28/05/2023

Regional Service Manager

www.DiaSys-Diagnostic.com

Contact: 022 3371 4323

Toll Free: 1800 120 1447

Registered Office

Plot No. A-821, T.T.C. Industrial Area,  
MIDC, Mahape, Navi Mumbai - 400710  
Maharashtra, India

**DiaSys Diagnostics India Private Limited**

The Specialists for diagnostic system solutions  
Point-of care testing and automated Laboratory systems

CIN : U74999MH2013PTC248371


Email : info@diasys.in

Tel. No. : +91-022-71068331


Toll Free : 1800 120 1447

www.diasys-diagnostics.com

### Light Quantity Checkup for all wavelength and found all in range

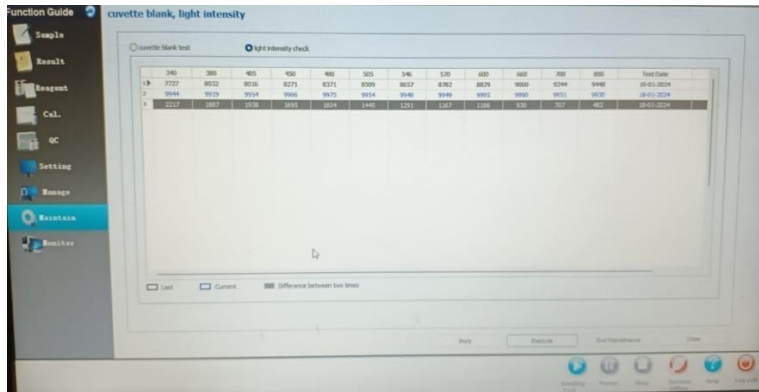

### Cuvette Blank test for all cuvettes found in range

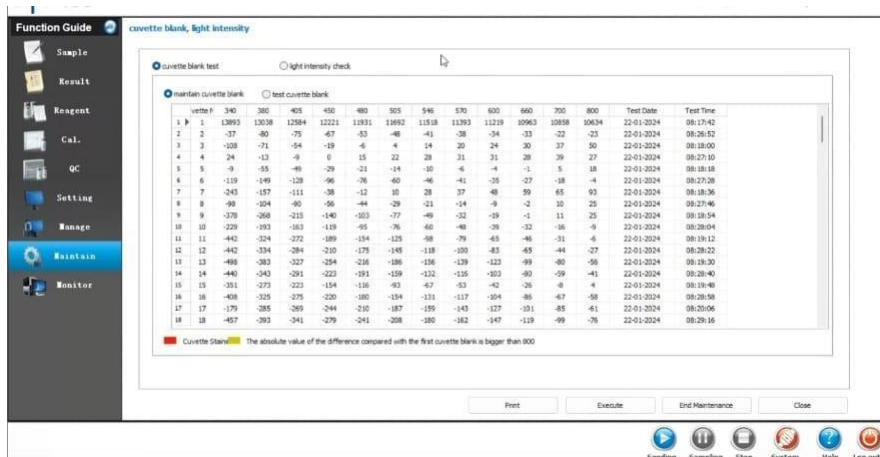


Asmat Sami

# Calibration Curve

Calibration Method: 2 Point Linear

Test Item: ALB


Calibrator: 1

Calibration Time: 25-07-2023 11:30:24

Absorbance: 0.1313

S1 Abs	K	A	B	C
0.1313	11.022	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.1298	6.	0	0
2.	3.64	0.45245	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

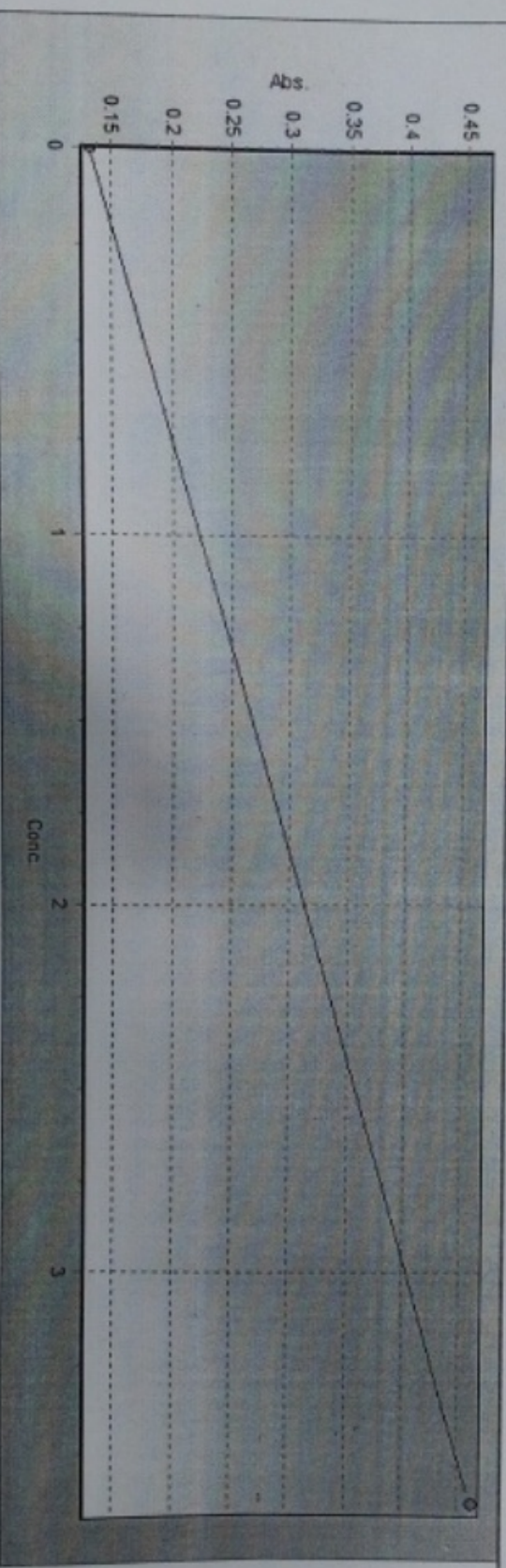

# Calibration Curve

Calibration Method: 2 Point Linear

Test Item: ALP

Calibrator: 1


Calibration Time: 27-05-2023 8:21:21

Absorbance: 9.99999999999998E

SI Abs	K	A	B	C
9.99999999999998E-5	5691.9643	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.0027	6.	0	0
2.	255	0.04615	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

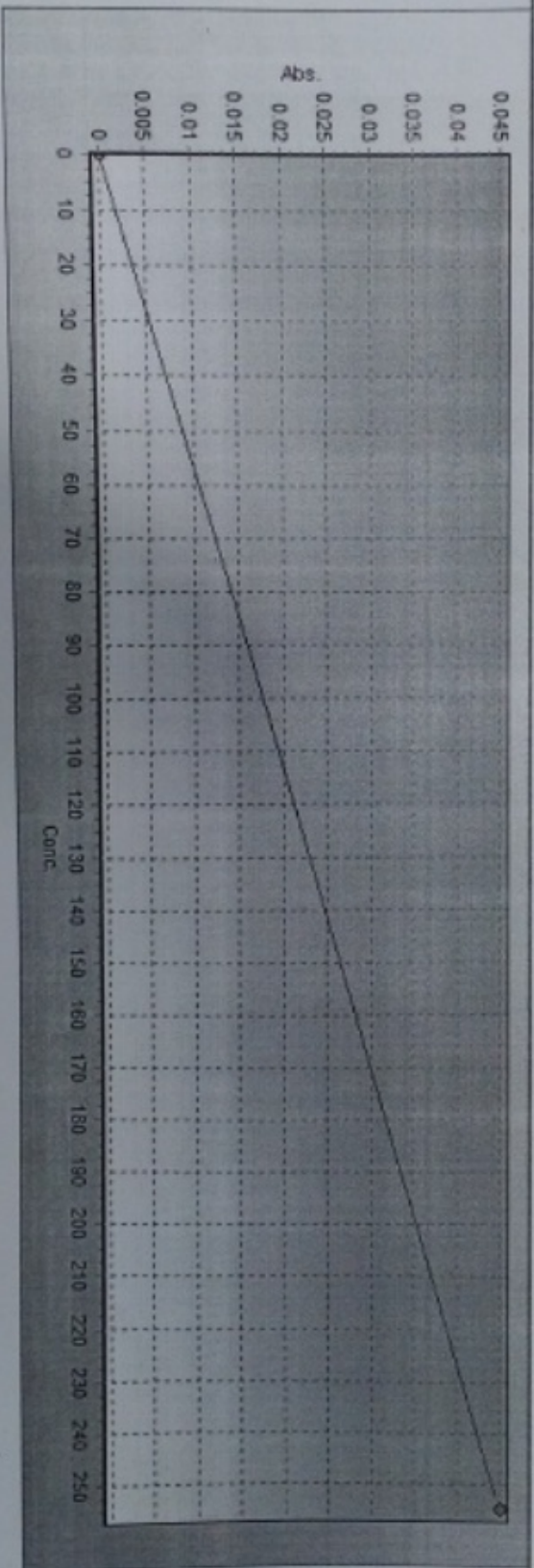

*[Handwritten signature]*

# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:21:52


Test Item: ALT

Absorbance: -0.0008

Calibrator: 1

	S1 Abs	K	A	B	C
	-0.0008	-3198.6644	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	-0.0011	6.	0	0
2.	95.8	-0.03275	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

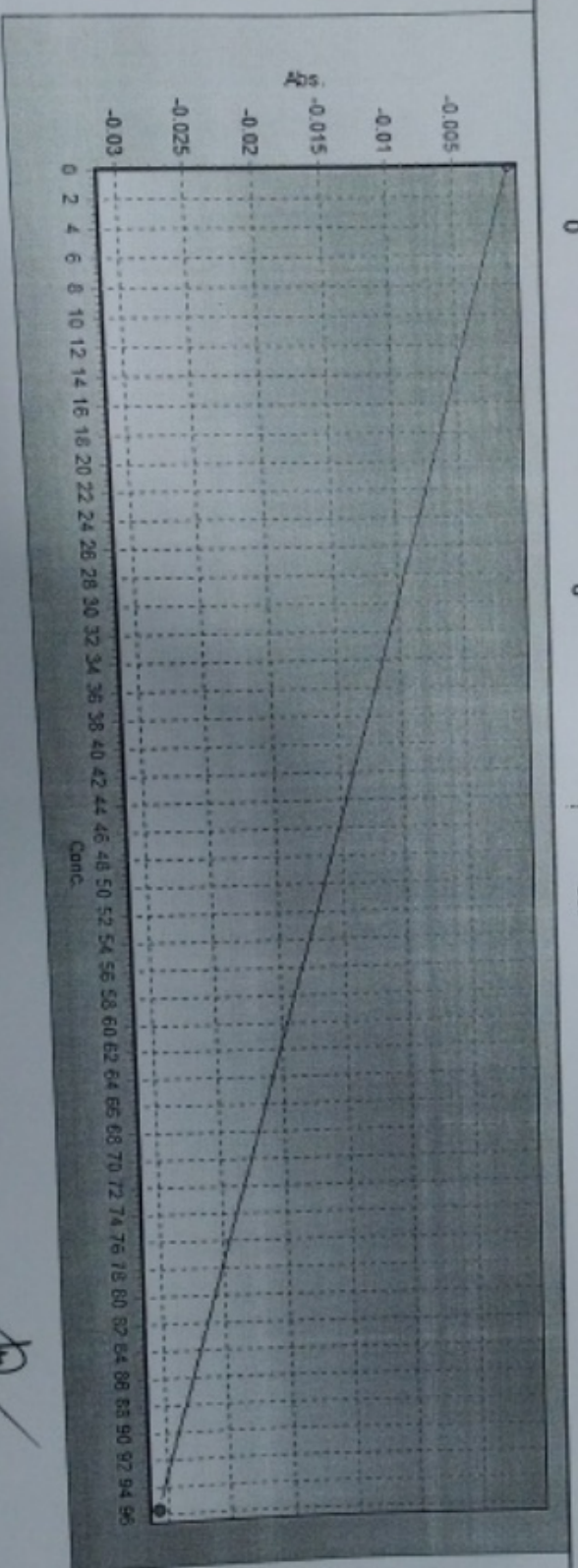

*Handwritten signature*

Print By: 1

Print Time: 05-03-2024 16:34:43

# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:22:52


Test Item: AST

Absorbance: -0.00165

Calibrator: 1

S1 Abs	K	A	B	C
-0.00165	-3267.9739	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	-0.00165	6.	0	0
2.	100	-0.03225	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

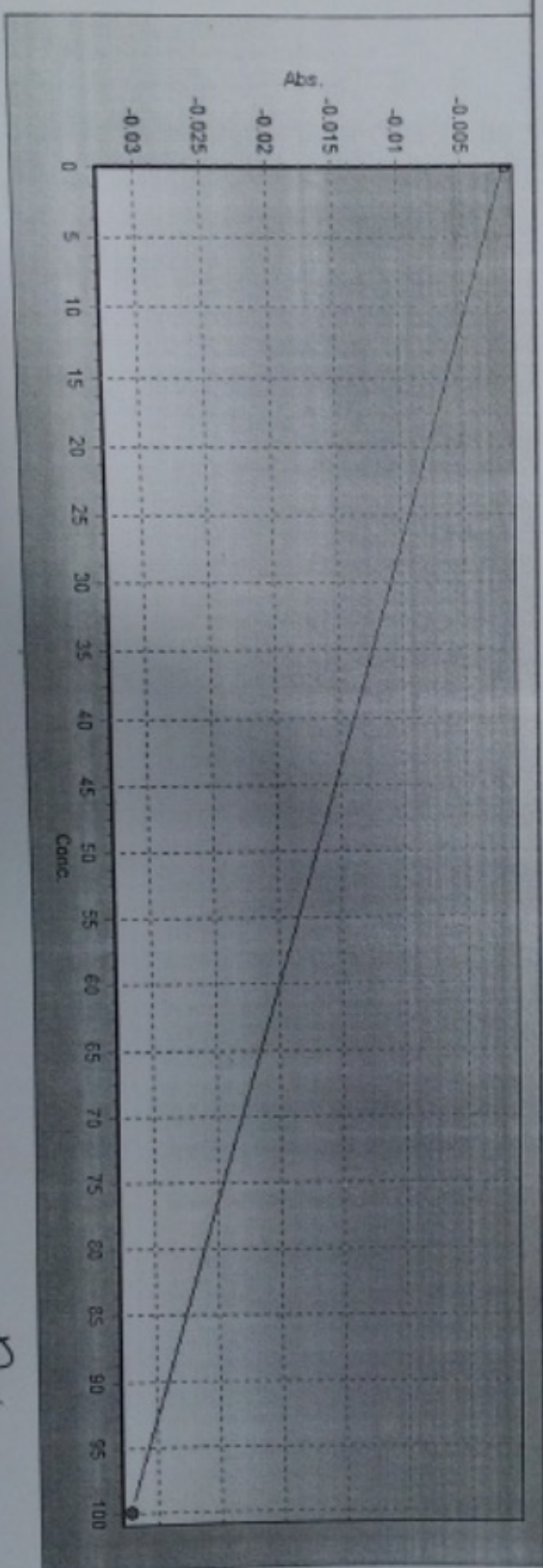

Print Bv: 1

Print Time: 05-03-2024 16:35:27

# Calibration Curve

Calibration Method: 2 Point Linear

Test Item: GREAT JAFFE


Calibrator: 1

Calibration Time: 27-05-2023 8:22:06

Absorbance: 0.00155

S1 Abs	K	A	B	C
0.00155	119.2488	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.00175	6.	0	0
2.	3.81	0.0375	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

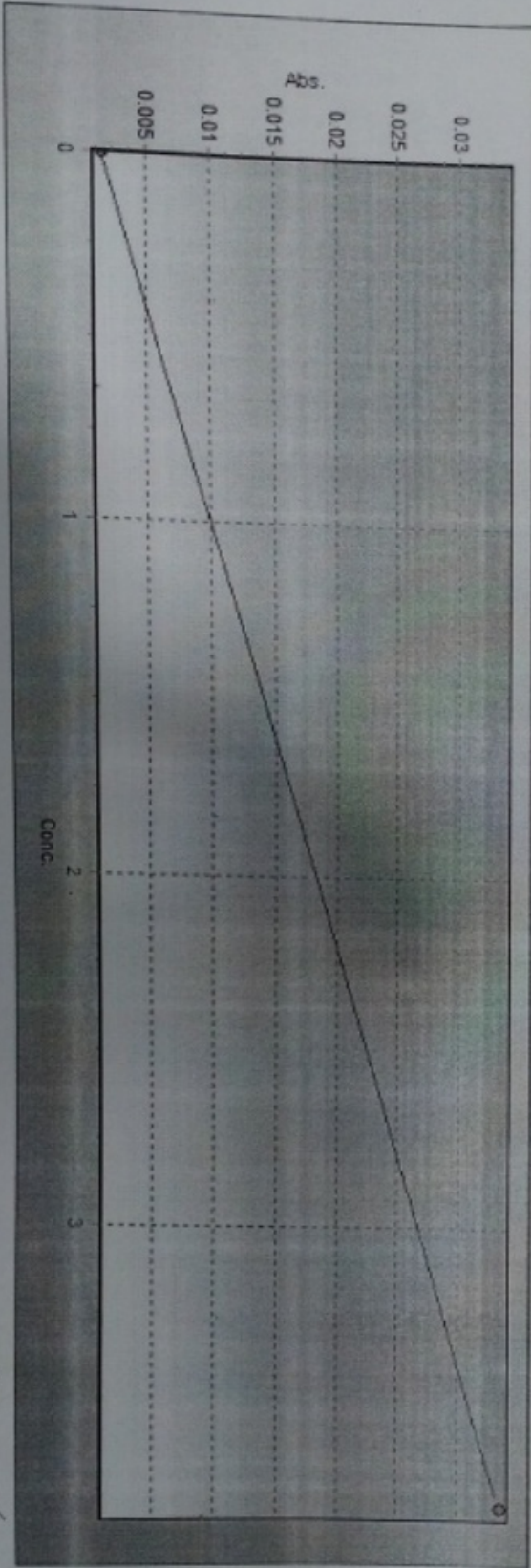

*[Handwritten Signature]*

# Calibration Curve


Calibration Method: 2 Point Linear  
 Test Item: DBIL  
 Calibrator: 1

Calibration Time: 27-05-2023 8:26:37  
 Absorbance: -5E

S1 Abs	K	A	B	C
-5E-5	41.7732	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	-5E-5	6.	0	0
2.	2.45	0.0586	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

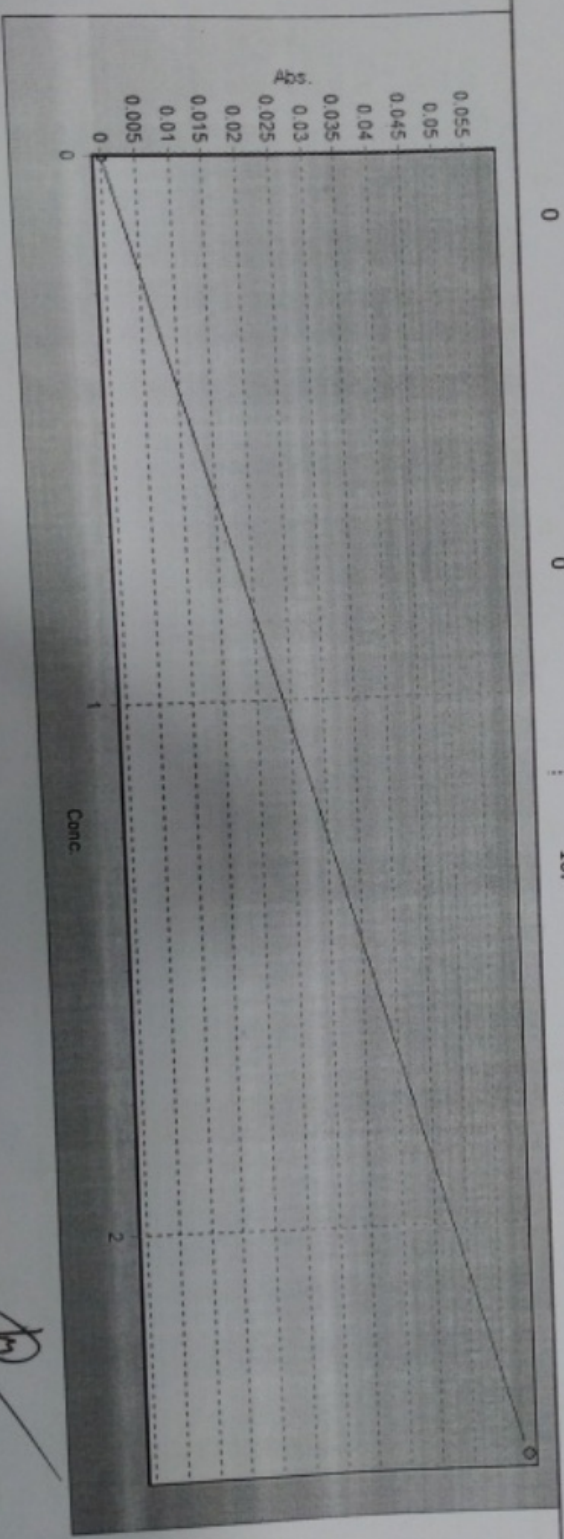

*[Handwritten Signature]*



# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:26:38


Test Item: GLU

Absorbance: 0.04045

Calibrator: 1

S1 Abs	K	A	B	C
0.04045	258.6915	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.04045	6.	0	0
2.	189	0.77105	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

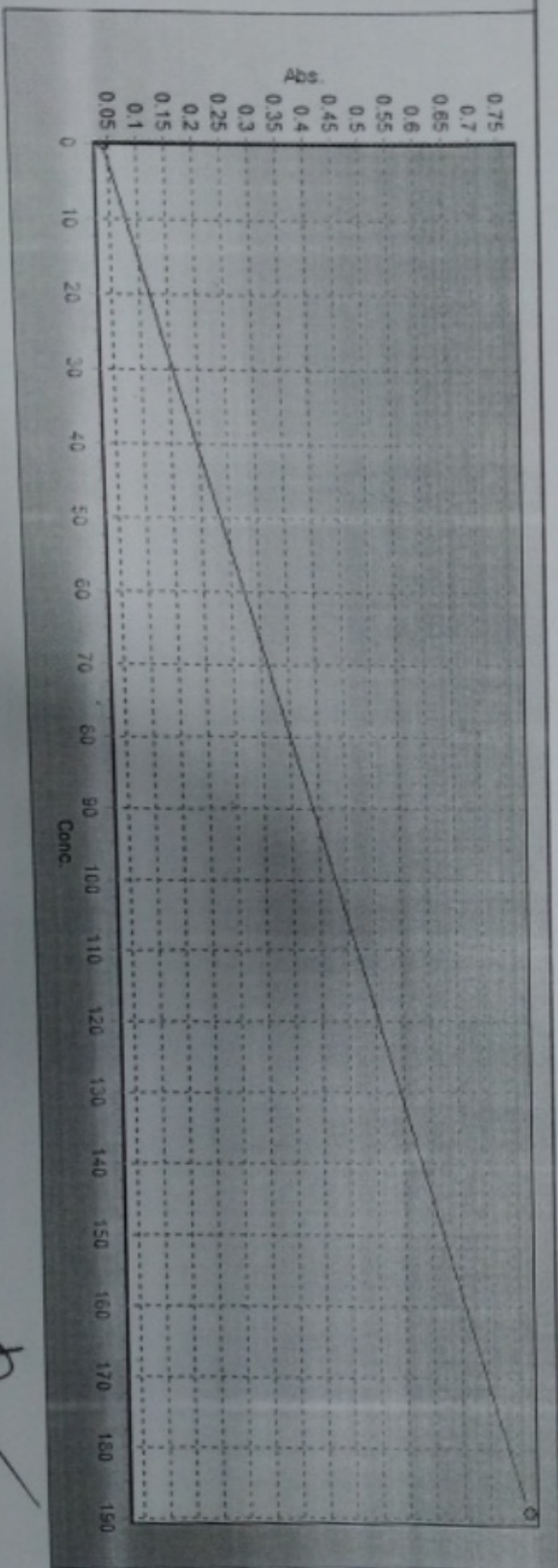

*[Handwritten signature]*

# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:28:37


Test Item: TBIL

Absorbance: 0.0088

Calibrator: 1

SI Abs	K	A	B	C
0.0088	88.6731	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.0088	6.	0	0
2.	4.11	0.05515	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

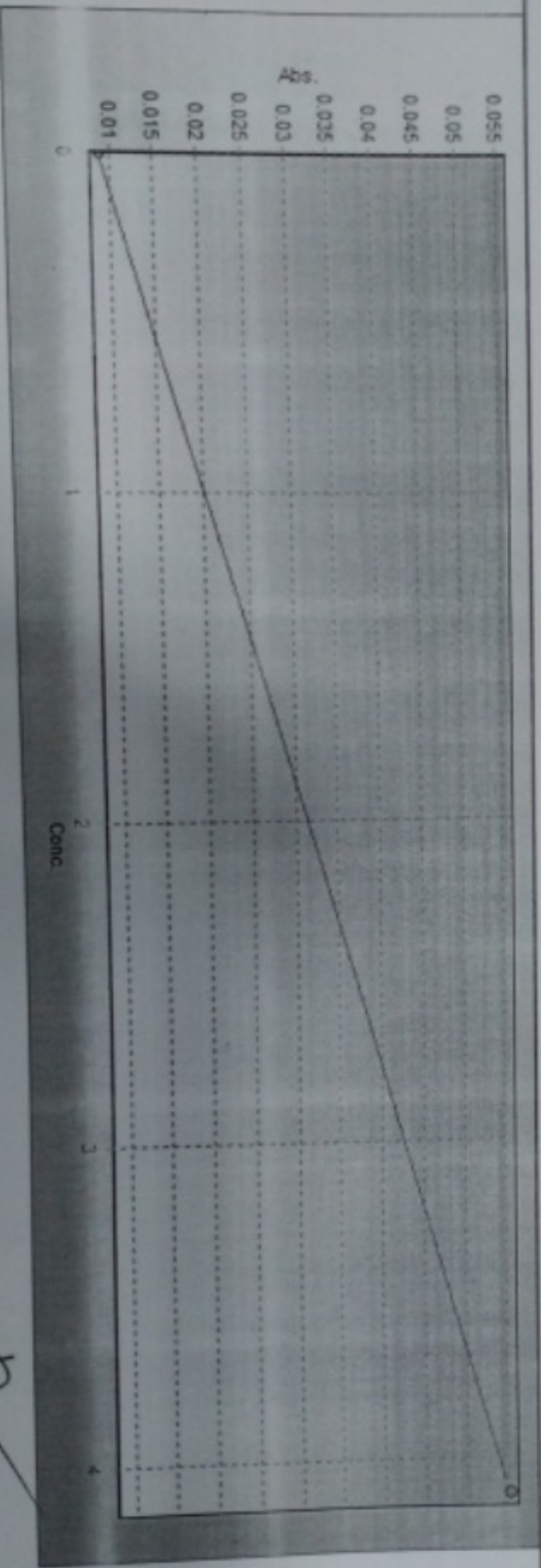

*[Handwritten Signature]*

# Calibration Curve

Calibration Method: 2 Point Linear

Test Item: TC


Calibrator: 1

Calibration Time: 27-05-2023 8:26:06

Absorbance: 0.02285

S1 Abs	K	A	B	C
0.02285	530.5919	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.0409	6.	0	0
2.	160	0.34115	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

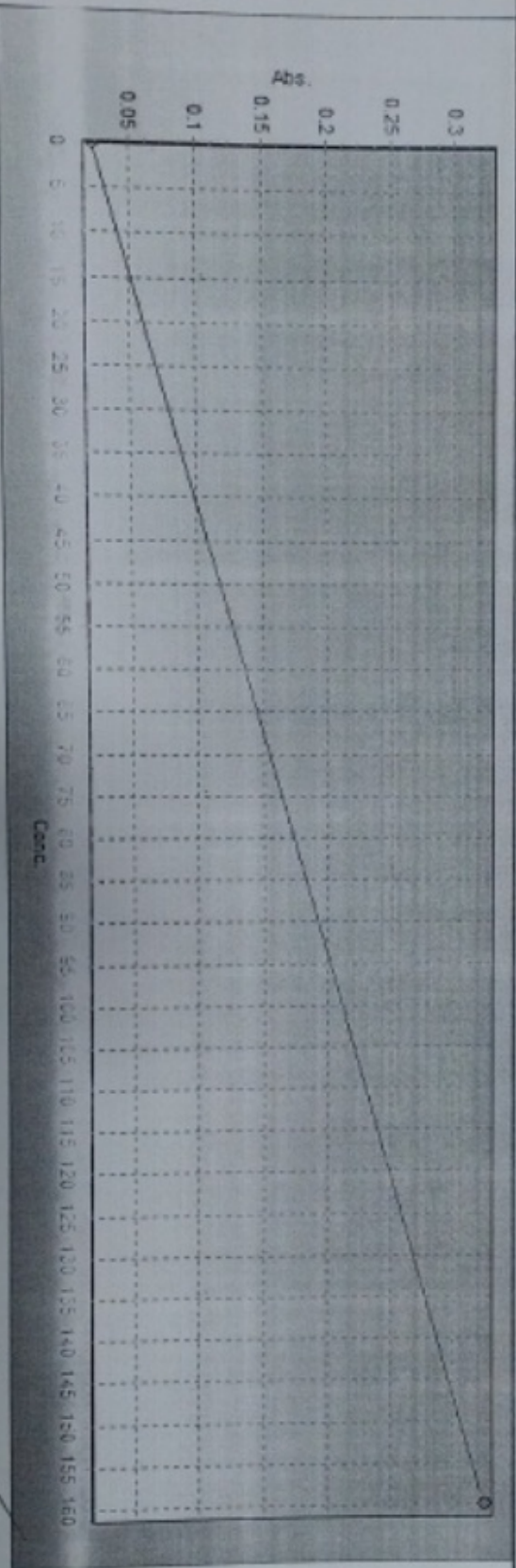

*(Handwritten signature)*

# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:27:06


Test Item: TG

Absorbance: 0.10925

Calibrator: 1

S1 Abs	K	A	B	C
0.10925	718.9216	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.10925	6.	0	0
2.	144	0.30955	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

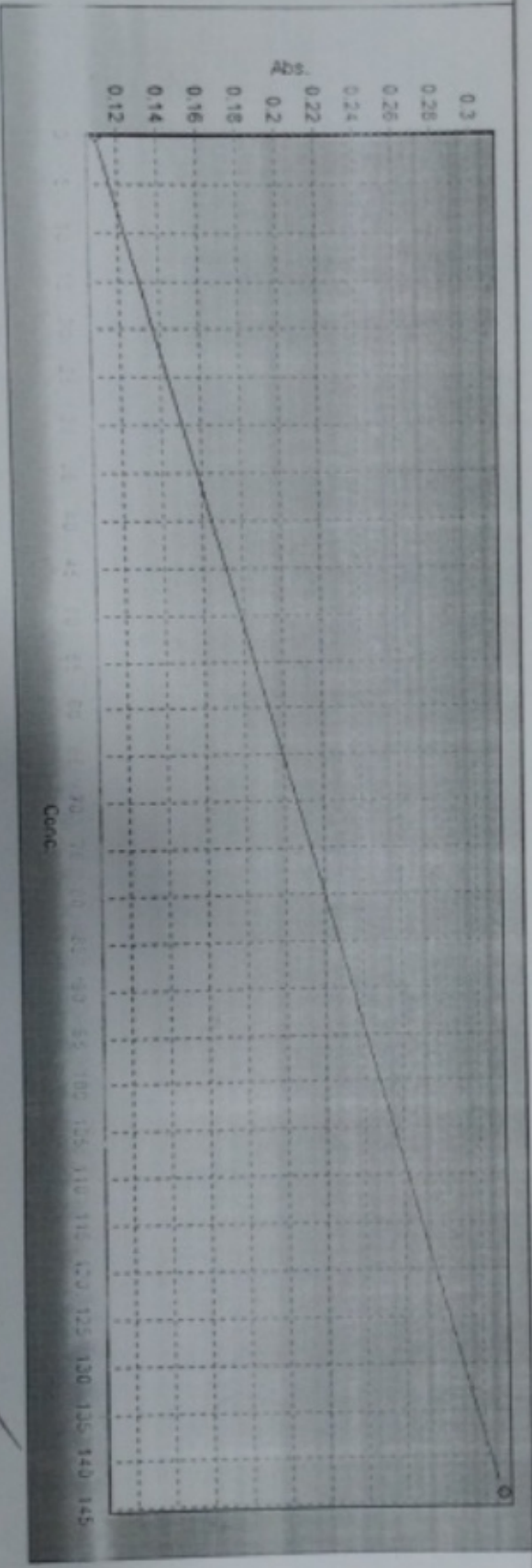

# Calibration Curve

Calibration Method: 2 Point Linear


Calibration Time: 27-05-2023 8:31:06

Test Item: TP

Absorbance: 0.00715

Calibrator: 1

	S1 Abs	K	A	B	C
	0.00715	43.2973	0	0	0
Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.00715	6.	0	0
2.	5.41	0.1321	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

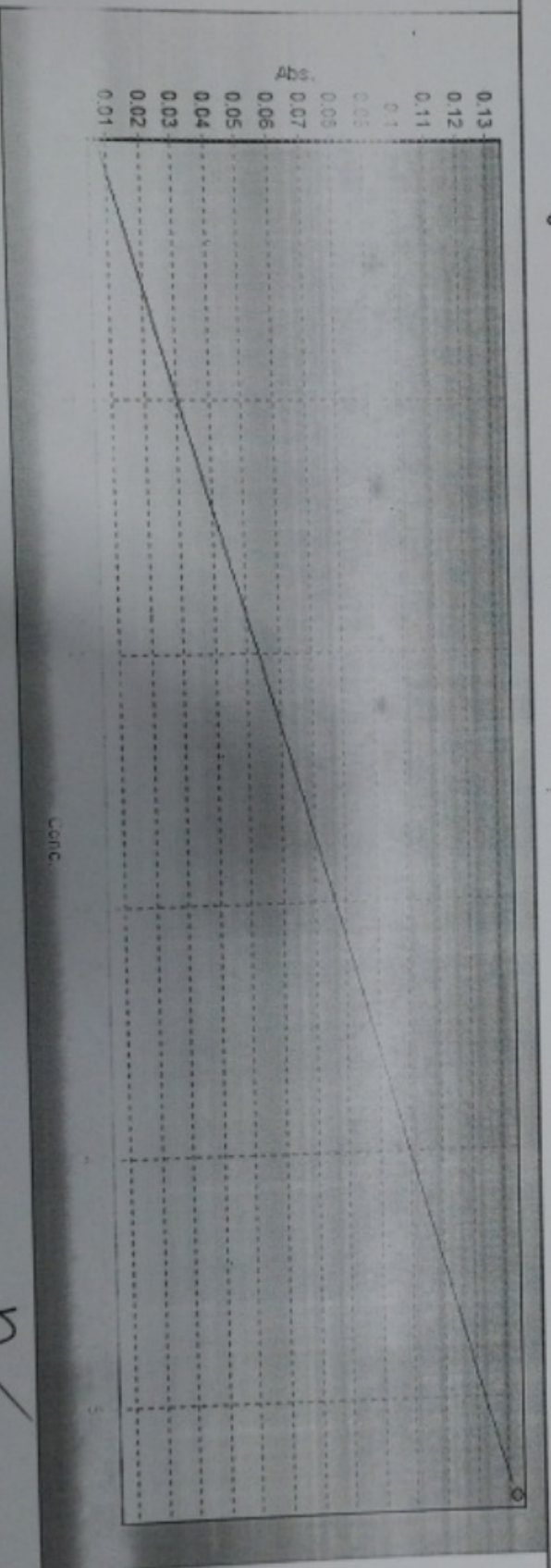

*[Handwritten Signature]*

# Calibration Curve

Calibration Method: 2 Point Linear

Test Item: UA


Calibrator: 1

Calibration Time: 27-05-2023 8:30:21

Absorbance: 0.003

	S1 Abs	K	A	B	C
	0.003	92.1902	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	0.003	6.	0	0
2.	5.43	0.0619	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

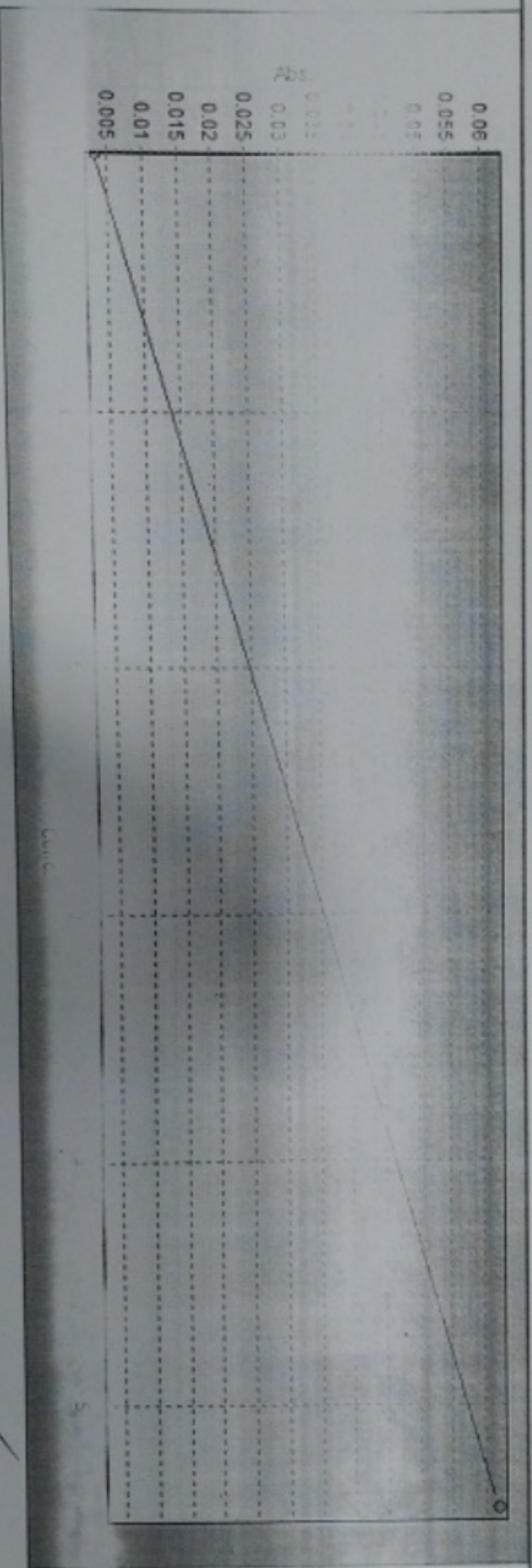

A

# Calibration Curve

Calibration Method: 2 Point Linear

Calibration Time: 27-05-2023 8:29:52


Test Item: UREA

Absorbance: -0.00965

Calibrator: 1

SI Abs	K	A	B	C
-0.00965	-1361.5928	0	0	0

Calibration	Concentration	Abs.	Calibration	Concentration	Abs.
1.	0	-0.0047	6.	0	0
2.	106	-0.0803	7.	0	0
3.	0	0	8.	0	0
4.	0	0	9.	0	0
5.	0	0	10.	0	0

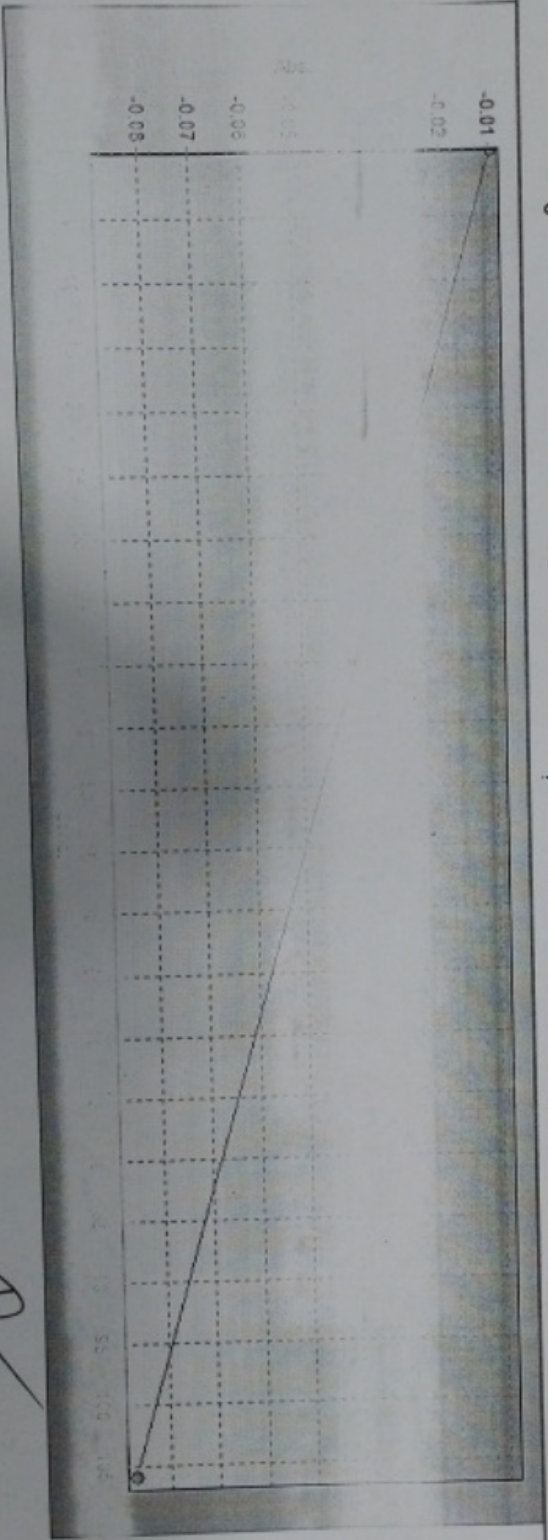

*[Handwritten Signature]*