



Installation Certificate for Standard F 200

This is to certify that Standard F 200 Serial No FA20FAITG1485 is successfully installed and commissioned at Lupin Diagnostics, Vijaywada . The installation protocol/checklist has been successfully completed for the above instrument.

For **SD Biosensor Healthcare Pvt. Ltd.**



Aamir Ayaz Khan

Assistant Product Manager

Installation Report**Customer Details****Instrument Details**

Name: Lupin Diagnostics
Contact Person(s) : Mr Prasad
Contact No : 9820355010
Address : Vijaywada

Model: Standard F 200
Serial Number : FA20FAITG1485
Installation Date : 03-08-2023

Accessories and Consumable Supplied

Power adapter: Yes
Calibration set: Yes
Test kit, Buffer: Yes
User instruction guide: Yes

Installation Checks


Environment: Closed
Power Adapter: Working
Temperature: 26 degree Celsius (in limits between 15 degree to 32 degree Celsius)
UPS: Yes

Training Imparted**Assays Performed**

Operating Procedure
Maintenance
Dos and Don'ts
Calibration

Customer Remarks**Installer's Remarks**

Customer's Signature

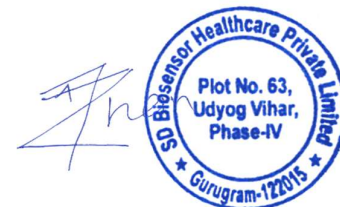


Installer's Signature

Operations Qualification for Standard F 200

1. Verified Cal 1 Calibration (blank check) : Done
2. Verified Cal 2 Calibration (UV light source and detector check) : Done
3. Verified Cal 3 Calibration (RBG light source and three level band Intensity (Strong positive, weak positive and negative) check) : Done

For **SD Biosensor Health Care Pvt. Ltd.**

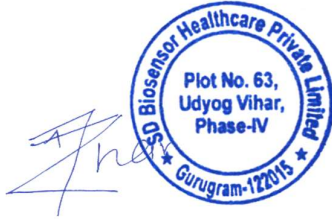

Aamir Ayaz Khan

Assistant Product Manager

Performance Certificate – Standard F 200



Instrument Certification:

Study data has determined that the instrument described in this document either meets all criteria outlined in this Performance Qualification Protocol or exceptional conditions have been identified and documentation included.

It is recommended to do Instrument calibration by calibration set once in 30 days for checking optical parts of analyzer. There is no preventive maintenance required in this equipment. SD Biosensor will update software of system free of cost as when required to run machine at optimum levels.

This instrument is ready for specific usage.

Protocol Performed By:



Aamir Ayaz Khan

Assistant Product Manager

In Field Calibration Report

Standard F200
Serial No : FA20FAITG1485

Operator
Name: S D Biosensor Healthcare Pvt Ltd

Instrument Owner
Name: Lupin Diagnostics, Vijaywada

Control Device	Lot No.	Expiry Date
CAL 1	6723B11F1BE02	31-10-2024
CAL 2	6723B11F1BE02	31-10-2024
CAL 3	6723B11F1BE02	31-10-2024
C -A1C – 1	9B24951AC	04-03-2025
C – A1C -2	9B24951AC	04-03-2025

Operating Condition	
Temperature	15-32 °C

Result of Control Device Measurements			
Control Device	Measured Value		Accepted Range
	1	2	
CAL 1	Ok	Ok	N/A
CAL 2	Ok	Ok	N/A
CAL 3	Ok	Ok	N/A
9E129V1AC – L1	5.4	5.3	4.3 % - 6.7 %
9E129V1AC – L2	9.1	9.2	7.5 % - 11.5 %

Final Result		
The Operator hereby confirm that the In Field Calibration was performed in accordance to the instructions given in " Standard F200 , In Field Calibration Procedure. This certificate is valid for 100 days.		
Date: 23-10-2024	Authorized Operator Signature	Instrument Responsible Signature

